

Pycnometry and Continuous Distillation of Ethanol-Water mixtures

Joseph Zimmerman

Abstract—In this paper, two experiments were conducted on ethanol-water mixtures. A pycnometer and quartz-tubed densitometer were used to determine liquid densities across a variety of compositions of ethanol and water. Alcohol-gauging tables from the US Alcohol and Tobacco Tax and Trade Bureau (TTB) were used to correlate theoretical densities of various compositions of ethanol and water at various temperatures to apparent proofs at 60F and finally correct to true proof and mole fraction of ethanol. Densities, measured through the pycnometer and densitometer, were compared.

Index Terms—Article submission, IEEE, IEEEtran, journal, L^AT_EX, paper, template, typesetting.

I. INTRODUCTION

PYCNOMETRY is a fundamental laboratory technique used to determine the density of solid or liquid substances by measuring their mass and volume and comparing them to a known volume. A pycnometer—a calibrated glass flask with a defined volume—is filled with the sample, and its mass is compared to the mass of the same flask filled with a reference fluid (water in our case). Digital densitometers can also be used to determine the density of a solution. Our digital densitometer measures density through a vibrating "U-Tube". The sample is introduced into a U-shaped borosilicate glass tube that is being excited to vibrate at its characteristic frequency electronically. The characteristic frequency changes depending on the density of the sample. Through determination of the characteristic

frequency the density of the sample can be calculated. In the first part of our experiment, we determined the densities of various compositions of ethanol-water through pycnometry and a digital densitometer. Alcohol-gauging tables from the US Alcohol and Tobacco Tax and Trade Bureau (TTB) were used to correlate theoretical densities of various compositions ranging from pure ethanol to pure water at various temperatures to apparent proofs at 60F and finally correct to extract true proof and mole fraction of ethanol.

This theoretical array of densities and corresponding mole fractions were then used in the second part of the lab, where we performed a continuous distillation of an ethanol-water mixture at two different reflux ratios, 4 and 5 in our case.

II. THEORY

A. Pycnometry

The pycnometer used in our experiment was a Kimble KiMax pycnometer. By measuring the mass of the empty pycnometer, the mass of the pycnometer filled with water (our reference fluid), the known density of water, and the mass of the pycnometer filled with the experimental solution, the density of the experimental solution can be determined by [1].

$$V = \frac{m_w - m_e}{\rho_{water}(t)} = \frac{m_{sample} - m_e}{\rho_{sample}} \quad (1)$$

where m_w is the mass of the pycnometer filled with water, m_e is the mass of the empty pycnometer, $\rho_{water}(t)$ is the density of water at temperature T, and ρ_{sample} is the density of the sample at T. This equation can be arranged to solve for ρ_{sample} .

$$\rho_{sample} = \frac{m_{sample} - m_{empty}}{m_{water} - m_{sample}} \cdot \rho_{water}(T) \quad (2)$$

B. Gauging alcohol with TTB tables

Table 6 from the TTB correlates specific gravities of ethanol-water mixtures at 60° F to proofs. In order to obtain proofs at other temperatures, one must convert the measured density (ρ_1) at a given temperature (T1) to the hypothetical density (ρ_2) that a glass hydrometer would indicate if calibrated at the standard 60°F reference temperature:

$$\rho_2 = \rho_1(1 + \alpha(T_1 - 60^\circ F)) \quad (3)$$

where α is the coefficient of thermal expansion for the hydrometer material (typically $25 \cdot \frac{10^{-6}}{^\circ C}$ for glass). At this corrected density, specific gravity can be calculated using the reference density of water at 60, °C 0.99904 g/cc. Then, Table 6 can be interpolated through a truncated Taylor series expansion to retrieve an apparent proof, C_2 , implemented as equation 1 in python:

$$C_2 = f(\gamma_1, T1) = [f + (\gamma_1 - \gamma_R) \frac{\partial f}{\partial \gamma}] \quad (4)$$

where γ_1 is the specific gravity of the sample corrected to 60°F. This apparent proof still has instrument and sample temperature effects, so it must be corrected to a true proof using TTB Table 1. The computed apparent proof C_2 is rounded down to the nearest integer and used as a reference for interpolation of Table 1.

$$C_3 = C_2 + (C_2 - C_R) \frac{\partial f}{\partial C} \Big|_{C_R, T_R} + (T_1 - T_R) \frac{\partial C}{\partial T} \Big|_{C_R, T_R} \quad (5)$$

where C_3 is the estimated true proof of the mixture corrected for temperature, T_1 . C_R and T_R are reference proofs and temperatures from values corresponding to the rounded-down C_2 in Table 1.

C. Error Propagation

1) *Pycnometer Masses*: Uncertainty for each of the 10 mixtures in the pycnometer, the pure water in the pycnometer, and the empty pycnometer were calculated on a 95% Confidence Interval (N=3).

2) *Pycnometer Densities*: Uncertainties for the mixture densities determined with the pycnometer propagates from uncertainties in mass measurements. The standard formula for error propagation was applied to equation (2)

$$\delta d_f = \sqrt{\left(\frac{\partial d_f}{\partial m_f} \delta m_f\right)^2 + \left(\frac{\partial d_f}{\partial m_e} \delta m_e\right)^2 + \left(\frac{\partial d_f}{\partial m_w} \delta m_w\right)^2 + \left(\frac{\partial d_f}{\partial T} \delta T\right)^2} \quad (6)$$

where δd_f is the uncertainty associated with the computed density, δm_f is the uncertainty associated with the sample-filled pycnometer mass measurements, δm_e is the uncertainty associated with the empty pycnometer mass measurements, and δm_w is the uncertainty associated with the water-filled pycnometer mass measurements and the partial with respect to temperature is estimated with finite differences

D. Distillation

1) *Vapor Liquid Equilibria*: The theoretical stages of our distillation column are assumed to be at vapor-liquid equilibrium (VLE). The VLE behavior of ethanol and water can be modeled using Antoine's equation for determining saturation pressures, and the Van Laar model for activity-coefficients. First, the activity coefficients, γ , are calculated

$$\gamma_1 = A_{12} \left(\frac{A_{21} X_2}{A_{12} X_1 + A_{21} X_2} \right)^2 \quad (7)$$

$$\gamma_2 = A_{21} \left(\frac{A_{12} X_1}{A_{12} X_1 + A_{21} X_2} \right)^2 \quad (8)$$

where A_{12} and A_{21} are Van Laar parameters for water and ethanol respectively. X_1 is the mole fraction of ethanol in solution and X_2 is the mole fraction of water. The Antoine Equation is used to determine saturation pressures of the pure components in the solutions and follows the form:

$$\log_{10} P^* = A - \frac{B}{T + C} \quad (9)$$

where P^* is the saturation pressure of the component and A, B, and C are tabulated constants. The VLE can then be modeled through the modified Raoult's Law:

$$\gamma_1 x P^{etOH*} = y P \quad (10)$$

$$\gamma_2 (1 - x) P^{H_2O*} = (1 - y) P \quad (11)$$

A solver like SciPy optimize's root function can be used to solve for saturation temperatures at a given pressure and a T_{xy} diagram can be constructed

2) *McCabe-Thiele diagram*: McCabe-Thiele analysis is based on the McCabe-Thiele assumption: that liquid is vaporized at the same rate vapor is condensed on a molar basis. With this assumption, two mole balances can be constructed: one for the operating (rectifying) line and one for the stripping line. These lines will be plotted against the T_{xy} equilibrium curve described earlier. The operating line passes through the point (x_p, x_p) , where x_p is the product ethanol mol fraction, and follows the equation:

$$y_i = \frac{O}{O+1} x_{i-1} + \frac{x_p}{O+1} \quad (12)$$

Where O is the reflux ratio of the column. The stripping line passes through the point (x_w, x_w) , where x_w is the waste ethanol mol fraction, and follows the equation:

$$y_i = \frac{O+F}{O+1} x_{i-1} + \frac{1-F}{O+1} x_w \quad (13)$$

Where F is the total molar flux per mole product.

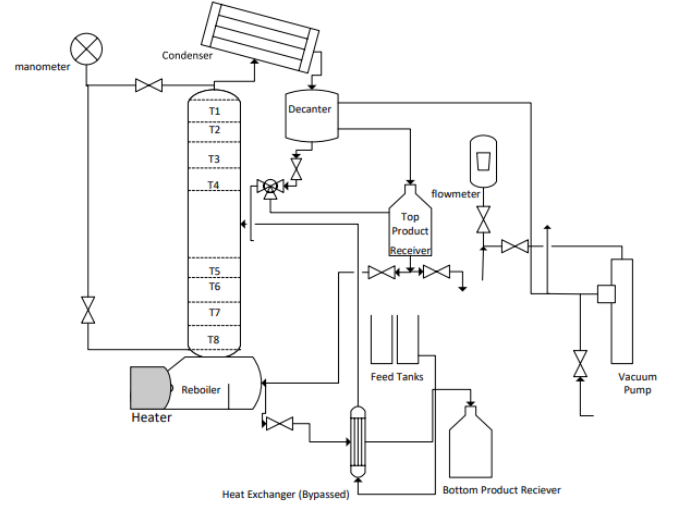


Fig. 1. Measured voltage readings from a PendoTec single-use pressure sensor recorded at various static pressures from a gas piston apparatus containing air. Static pressures were induced by adding weights of $[84.248 \pm 0.007\text{g}; 164.817 \pm 0.007\text{g}; 239.434 \pm 0.002\text{g}; 328.00 \pm 0.01\text{g}$ and $402.531 \pm 0.003\text{g}]$ to the apparatus and appear in the figure from left to right. Error bars were calculated based on a 95% confidence interval of recorded voltage readings in mV across three trials for each pressure (N=2). Uncertainty

III. EXPERIMENTAL METHODS

A. Pycnometry

1) *Pycnometer Calibration*: The pycnometer was rinsed thoroughly with ethanol and dried with compressed air. The mass of the empty (stopped and capped) pycnometer was weighed 3 times with an analytical balance and recorded as M_e . The pycnometer was emptied, cleaned following the above procedure, and filled with deionized water (DIW). The mass of the water-filled pycnometer was weighed 3 times with an analytical balance and recorded as M_w . Then, 3 readings of density and temperature were taken with the densitometer and recorded as ρ_w .

2) *Pycnometer Sample Measurements*: The pycnometer was emptied, rinsed thoroughly with ethanol and dried with compressed air. 10 ethanol-water mixtures ranging from pure ethanol to pure water were constructed. Various weights of ethanol and water were measured over 3 trials then combined in a beaker to form the samples we measured, cleaning the beaker and pycnometer before switching concentrations. The mass of the empty (stopped and capped) pycnometer was

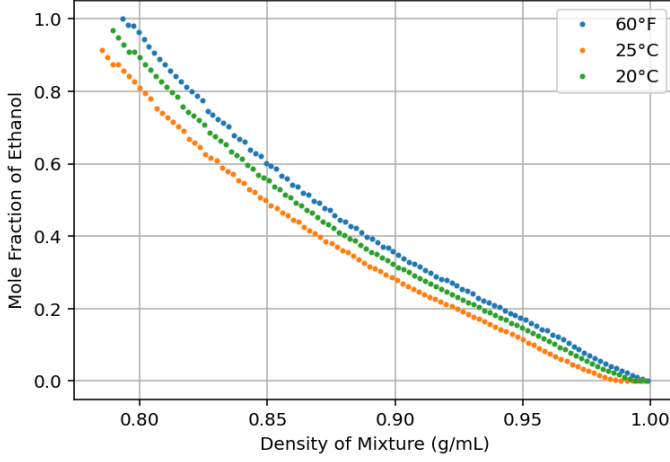


Fig. 2. Measured voltage readings from a PendoTec single-use pressure sensor recorded at various static pressures from a gas piston apparatus containing air. Static pressures were induced by adding weights of $[84.248 \pm 0.007\text{g}; 164.817 \pm 0.007\text{g}; 239.434 \pm 0.002\text{g}; 328.00 \pm 0.01\text{g}$ and $402.531 \pm 0.003\text{g}]$ to the apparatus and appear in the figure from left to right. Error bars were calculated based on a 95% confidence interval of recorded voltage readings in mV across three trials for each pressure (N=2). Uncertainty

weighed 3 times with an analytical balance and recorded as M_e . Then, 3 readings of density and temperature were taken at every concentration with the densitometer and recorded as ρ_{dN} .

B. Continuous Distillation

IV. RESULTS

A. Calibration

The calibration value was determined from the slope of the linear regression of recorded voltages, in mV, against varied pressures, in psi. This value was divided by 9V to match PendoTec specifications to be $0.270756 \pm 0.0000001 \frac{\text{mV}}{\text{psi} \cdot 9\text{V}}$.

This value was used to determine the pressure in the apparatus from voltage readings in the following Ruchardt experiments

Note that often IEEE papers with multi-part figures do not place the labels within the image itself (using the optional argument to `\subfloat[]`), but instead will reference/describe all of them (a), (b), etc., within the main caption. Be aware that for `subfig.sty` to generate the (a), (b), etc., subfigure labels, the optional argument to `\subfloat` must be present.

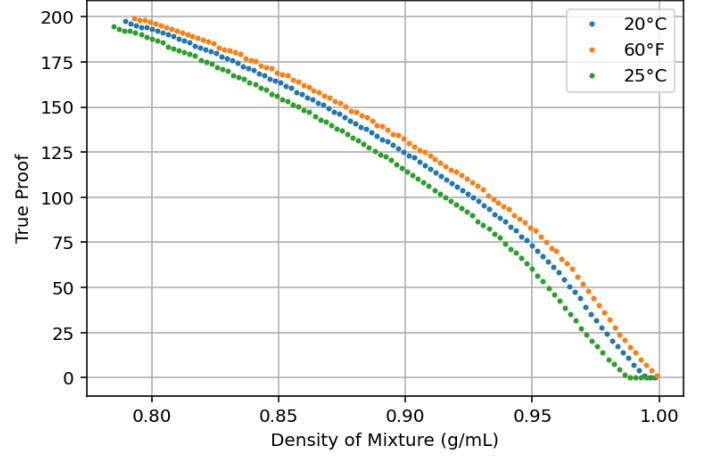


Fig. 3. Measured voltage readings from a PendoTec single-use pressure sensor recorded at various static pressures from a gas piston apparatus containing air. Static pressures were induced by adding weights of $[84.248 \pm 0.007\text{g}; 164.817 \pm 0.007\text{g}; 239.434 \pm 0.002\text{g}; 328.00 \pm 0.01\text{g}$ and $402.531 \pm 0.003\text{g}]$ to the apparatus and appear in the figure from left to right. Error bars were calculated based on a 95% confidence interval of recorded voltage readings in mV across three trials for each pressure (N=2). Uncertainty

If a subcaption is not desired, leave its contents blank, e.g., `\subfloat[]`.

V. DISCUSSION

Note that, for IEEE-style tables, the `\caption` command should come BEFORE the table. Table captions use title case. Articles (a, an, the), coordinating conjunctions (and, but, for, or, nor), and most short prepositions are lowercase unless they are the first or last word. Table text will default to `\footnotesize` as the IEEE normally uses this smaller font for tables. The `\label` must come after `\caption` as always.


VI. TECHNICAL ADVANCES

VII. REFERENCES SECTION

You can use a bibliography generated by BibTeX as a .bbl file. BibTeX documentation can be easily obtained at: <http://mirror.ctan.org/biblio/bibtex/contrib/doc/> The IEEEtran BibTeX style support page is: <http://www.michaelshell.org/tex/ieeetran/bibtex/>

TABLE I
GAS PROPERTIES

Gas	Condition	P_0 (Pa)	C (Pa)	κ (s ⁻¹)	ω (rad s ⁻¹)	ϕ (rad)
Air	Empty	102171.80 \pm 7.20	4518.24 \pm 886.00	37.53 \pm 2.97	217.85 \pm 2.87	-9.62 \pm 3.02
	1 Weight	103076.35 \pm 2.59	2061.26 \pm 545.07	14.59 \pm 1.03	118.06 \pm 0.64	-10.88 \pm 7.99
	2 Weight	103389.06 \pm 1.81	4264.83 \pm 513.42	10.45 \pm 0.76	90.98 \pm 0.63	-8.53 \pm 6.94
	3 Weight	104821.81 \pm 12.52	-1212.16 \pm 513.84	10.03 \pm 0.56	79.06 \pm 0.73	-15.93 \pm 9.20
Argon	Empty	102269.40 \pm 4.39	-1487.97 \pm 3403.32	16.59 \pm 4.37	207.72 \pm 0.38	-12.77 \pm 6.63
	1 Weight	103093.24 \pm 3.53	23.81 \pm 945.56	15.69 \pm 0.08	120.07 \pm 0.51	-15.62 \pm 4.69
	2 Weight	104000.11 \pm 21.84	1201.48 \pm 2621.70	17.28 \pm 1.22	80.95 \pm 0.35	-31.05 \pm 8.21
	3 Weight	104882.66 \pm 36.64	349.51 \pm 7984.72	14.13 \pm 0.50	77.39 \pm 0.91	-15.90 \pm 11.16
CO ₂	Empty	102304.81 \pm 15.52	882.70 \pm 1867.09	32.71 \pm 5.98	178.53 \pm 1.66	-22.61 \pm 15.44
	1 Weight	102216.12 \pm 3.98	890.54 \pm 3907.62	14.78 \pm 3.93	106.12 \pm 2.15	-22.41 \pm 23.91
	2 Weight	103976.97 \pm 3.03	-1087.62 \pm 2016.17	9.09 \pm 1.17	83.14 \pm 0.42	-25.64 \pm 13.47
	3 Weight	104000.26 \pm 14.20	-6199.73 \pm 518.83	6.59 \pm 1.27	68.67 \pm 0.14	-13.64 \pm 11.07

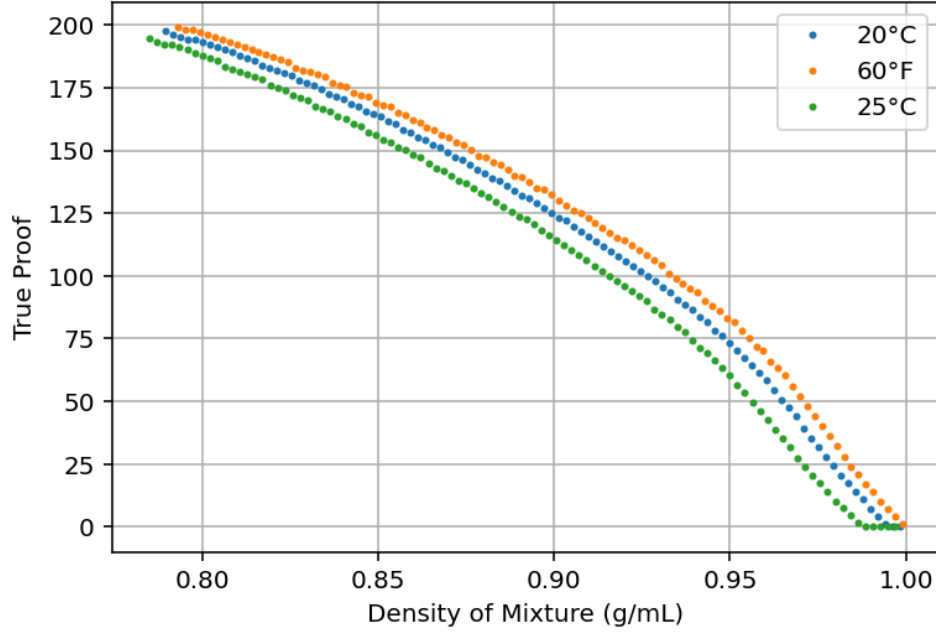

(a)

Fig. 4. Dae. Ad quatur autat ut porepel itemoles dolor autem fuga. Bus quia con nessunti as remo di quatus non perum que nimus. (a) Case I. (b) Case II.

TABLE II
 κ VALUES FOR AIR, ARGON, AND CO₂

Gas	κ_{Empty}	$\kappa_{1 \text{ Weight}}$	$\kappa_{2 \text{ Weight}}$	$\kappa_{3 \text{ Weight}}$
Air	1.37 \pm 4 $\times 10^{-21}$	0.95 \pm 3 $\times 10^{-21}$	1.10 \pm 3 $\times 10^{-21}$	1.20 \pm 3 $\times 10^{-21}$
Argon	1.74 \pm 3 $\times 10^{-21}$	1.61 \pm 4 $\times 10^{-21}$	1.53 \pm 2 $\times 10^{-21}$	1.64 \pm 3 $\times 10^{-21}$
CO ₂	1.28 \pm 1 $\times 10^{-21}$	1.36 \pm 4 $\times 10^{-21}$	1.32 \pm 2 $\times 10^{-21}$	1.28 \pm 2 $\times 10^{-21}$

VIII. SIMPLE REFERENCES

REFERENCES

You can manually copy in the resultant .bbl file and set second argument of `\begin` to the number of references (used to reserve space for the reference number labels box).


- [1] *Mathematics Into Type*. American Mathematical Society. [Online]. Available: <https://www.ams.org/arc/styleguide/mit-2.pdf>
- [2] T. W. Chaundy, P. R. Barrett and C. Batey, *The Printing of Mathematics*. London, U.K., Oxford Univ. Press, 1954.

- [3] F. Mittelbach and M. Goossens, *The L^AT_EX Companion*, 2nd ed. Boston, MA, USA: Pearson, 2004.
- [4] G. Grätzer, *More Math Into L^AT_EX*, New York, NY, USA: Springer, 2007.
- [5] M. Letourneau and J. W. Sharp, *AMS-StyleGuide-online.pdf*, American Mathematical Society, Providence, RI, USA, [Online]. Available: <http://www.ams.org/arc/styleguide/index.html>
- [6] H. Sira-Ramirez, “On the sliding mode control of nonlinear systems,” *Syst. Control Lett.*, vol. 19, pp. 303–312, 1992.
- [7] A. Levant, “Exact differentiation of signals with unbounded higher derivatives,” in *Proc. 45th IEEE Conf. Decis. Control*, San Diego, CA, USA, 2006, pp. 5585–5590. DOI: 10.1109/CDC.2006.377165.
- [8] M. Fliess, C. Join, and H. Sira-Ramirez, “Non-linear estimation is easy,” *Int. J. Model., Ident. Control*, vol. 4, no. 1, pp. 12–27, 2008.
- [9] R. Ortega, A. Astolfi, G. Bastin, and H. Rodriguez, “Stabilization of food-chain systems using a port-controlled Hamiltonian description,” in *Proc. Amer. Control Conf.*, Chicago, IL, USA, 2000, pp. 2245–2249.

IX. BIOGRAPHY SECTION

If you have an EPS/PDF photo (graphicx package needed), extra braces are needed around the contents of the optional argument to biography to prevent the LaTeX parser from getting confused when it sees the complicated `\includegraphics` command within an optional argument. (You can create your own custom macro containing the `\includegraphics` command to make things simpler here.)

If you include a photo:

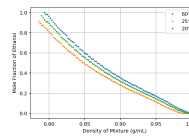

Michael Shell Use `\begin{IEEEbiography}` and then for the 1st argument use `\includegraphics` to declare and link the author photo. Use the author name as the 3rd argument followed by the biography text.

If you will not include a photo:

John Doe Use `\begin{IEEEbiographynophoto}` and the author name as the argument followed by the biography text.